


Phizzle

Digital Air Series

The only software-defined solution to automate data capture and remotely operate air particle counters

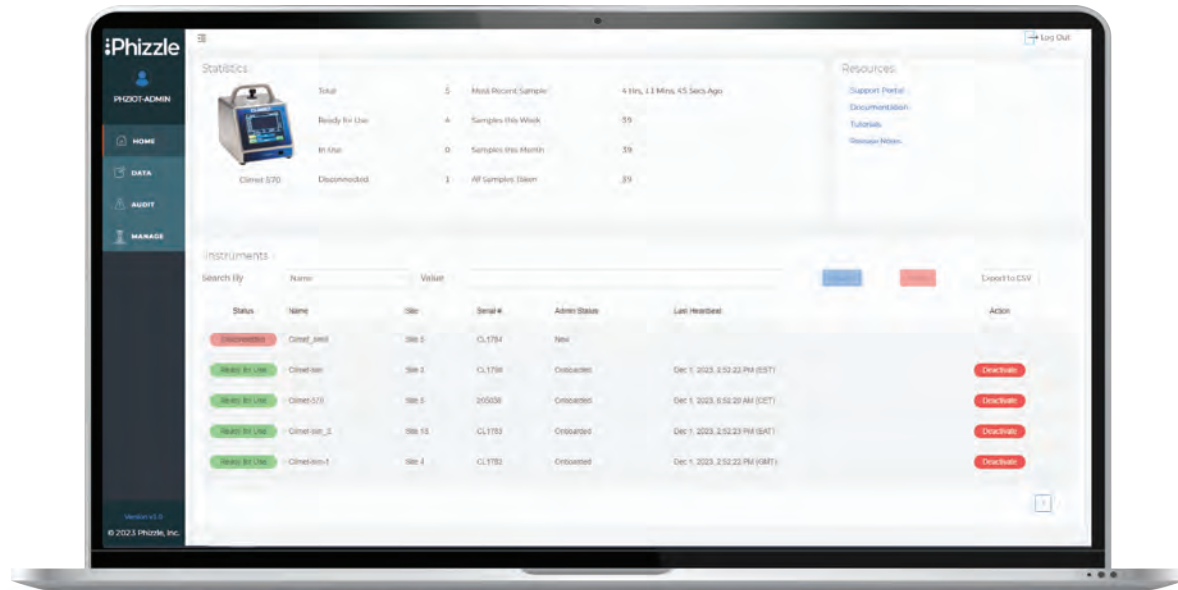

First to Market

Phizzle's channel partners are industry leaders in network infrastructure and instrument qualification. Together, we are building the Digital Air ecosystem to bring modern software benefits to cleanroom environments — benefits like digital record keeping, e-signature, remote operations, and automation. No other partnership in the industry can offer this type of comprehensive digital transformation.

Innovation: First of its kind

Digital Air is the only software-defined solution to automate data capture and remotely operate air particle counters used GMP/ISO applications.

Value: Guaranteed ROI

Any GMP/ISO customer with multiple air particle counters in multiple locations is guaranteed operational cost savings.

Proof of Value: In Production

The same software: FDA-approved, in-production with one of the most important drugs in the world.

Contact us for a demo:



www.phizzle.com

www.phizzle.com/request-a-demo

Phizzle

APC

Brands We Work With:



Climet



Beckman Coulter



Lasair



AeroTrak

It's Easy

Digital Air is so easy to use, we asked our intern to show you how it's done



Scan to watch the unboxing & setup video



It's Battle Tested!

The Digital Air platform was built from Phizzle's years of experience in the pharmaceutical manufacturing space. Presently, this technology is deployed at a Fortune 100 Pharma company making several of the world's most important drugs. This is an FDA-approved, GMP-ready solution that has passed the highest-level audits in pharmaceutical manufacturing. Compliance, security, and traceability are all in the product's DNA. Phizzle already operates in mission-critical environments across the globe and can be trusted to integrate with your processes and requirements to whatever use case suits you.

Phizzle Quick Start Guide Climet CI-x70

Welcome to Phizzle's Test Drive Program. You've received a box from our Sales team and are ready to begin using the Digital Air Series Platform with your Climet CI-x70. This Quick Start Guide will teach you how to onboard air particle counters to the platform — the first step to a fully digitized air particle counter workflow.

- Power on your Climet. Log into the Climet and enter the Main Menu.
- Swipe down from the Main Menu to enter the Settings page. Enter the WiFi Menu labeled RH in the upper left.
- Click to Enable WiFi.
- Click to see you have successfully connected — you should have a green WiFi signal in the upper left.
- Add the SSID of your WiFi network. Apply your Security settings. Reboot the Climet.
- Check to see you have successfully connected — you should have a green WiFi signal in the upper left.
- Enter the Main Menu again and click the Security icon, shaped like a key.
- Add an Admin user with both Name and Pass enabled. Make note of the password you set — this should be the same for all your air particle counters.
- Enter the Main Menu again and enter General Settings.
- Make note of the WiFi IP. Each air particle counter will be different; in this example, it is "192.168.8.158".

Scan to download the Quick Start Guide



www.phizzle.com